


REMODELLING
FACTOR



SLOUGHY/
EXUDATING
WOUNDS

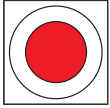

CURASORB*

Calcium-sodium alginate dressing

Optimises
fluid management



tyco
Healthcare

KENDALL



Advanced
Wound Care™

CURASORB*

Composition and Characteristics

Soft, sterile, non-woven dressing composed of calcium-sodium alginate fibres. In contact with wound exudate, develops the characteristic gel-like consistency of calcium alginate dressings. CURASORB remains intact in the presence of excessive wound fluid, and provides a gel block feature that minimises the spreading of the moisture to the outer perimeter of the dressing.

Product Features

Highly Absorbent	Absorbs approximately 15-20 times its weight in exudate Ideal for medium to heavily draining wounds Reduces dressing changes
Vertical Wicking	Low Sodium Content Enhances moist wound healing Binds exudate thus preventing skin maceration
Stable Gel Formula	Enhanced dressing integrity due to a higher concentration of Gluronic Acid (68%) than Mannuronic Acid (32%) Easy to locate dressing edges
Non-Adherent	Easy, painless dressing removal, minimises tissue trauma and maximises patient comfort
Conformable	Allows full mobility of dressing site
High Calcium Content	Aids haemostasis
Dressing and Rope Format	Can be customised / trimmed to fit a variety of wounds and wound sizes Ability to fill deep or tunnel wounds
CURASORB PLUS	Enhanced absorptive capabilities



Decubitus ulcer of the sacrum, 6 months old



Dressing with Curasorb



Result after 4 weeks

Indications

- Granulation tissue
- Moderately to heavy draining wounds
- Deep wounds with cavities
- Bleeding wounds
- Dehisced surgical wounds
- Venous ulcers
- Diabetic ulcers
- Arterial ulcers
- Donor sites
- Second degree burns
- Abrasions / lacerations / skin tears
- Pressure ulcers
- Other external wounds with moderate to heavy exudate

Ordering Information

CURASORB*

	REF	Pack Quantity	Case Quantity
30cm Rope	9231	5	20
5cmx5cm	9232	10	100
10cmx10cm	9233	10	50
10cmx10cm Plus	9236	10	50
20cmx10cm	9238	5	50
15cmx25cm	9239	10	50
10cmx14cm	9240	10	50
30cmx61cm	9242	5	5
61cm Rope	9243	5	20
91cm Rope	9244	5	20

Test data and study details on file. Reproduced with kind permission of Prof. Elia Ricci, University of Turin, ITALY
* CURASORB is a trademark of Beiersdorf AG. © 2003 Tyco Healthcare Group LP. All rights reserved.

P-WC-CURASO/GB

tyco
Healthcare

KENDALL

Tyco Healthcare UK Limited
154 Fareham Road
GOSPORT
Hampshire
PO13 0AS
UK
Tel: +44 (0) 1329 224226
Fax: +44 (0) 1329 224334

Tyco Healthcare Nederland BV
Hogeweg 105
530 1 LL ZALTBOMMEL
The NETHERLANDS
Tel: +31 418 57 66 00
Fax: +31 418 57 67 93
TycoHealthNed@TycoHealth.com
www.tycohealthcare.nl

Tyco Healthcare Regional Export
154 Fareham Road
GOSPORT
Hampshire
PO13 0AS
UK
Tel: +44 (0) 1329 224258
Fax: +44 (0) 1329 224086

Tyco Healthcare Norden AB
Gårdsvägen 18
SE-169 67 SOLNA
SWEDEN
Tel: +46 (0) 8 585 605 05
Fax: +46 (0) 8 585 605 06

Tyco Healthcare
Birchwood Court
Montrose Road
P.O. Box 8108
Halfway House
1685 MIDRAND
SOUTH AFRICA
Tel: +27 (0)11 265 2000
Fax: +27 (0)11 265 2050

www.tycohealthcare.com